




biodirect.

# sequência. fresagem

Velocidade de fresagem [RPM]	800 1000	800 1000	800 1000	800 1000	800 1000	800 1000
<b>Fresas</b>	Lança #2,2 <b>5088</b>	Helicoidal #2,2 <b>5073</b>	Cônica #2,8 <b>5084</b>	Cônica #3,2 <b>5085</b>	Cônica #3,6 <b>5096</b>	Cônica #4,0 <b>5116</b>

Osso tipo I e II

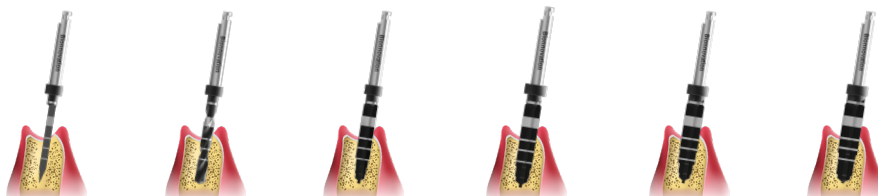

**Implantes**  
Diâmetro

3,75  
4,00

3,75  
4,00

3,75  
4,00

3,75  
4,00

3,75  
4,00

4,00  
Opcional

Osso tipo III e IV



**Implantes**  
Diâmetro


3,75  
4,00

3,75  
4,00

3,75  
4,00

3,75  
4,00

4,00

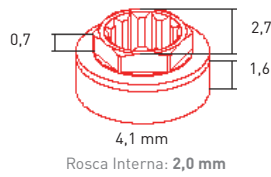


biodirect.3.3

# sequência. fresagem

Velocidade de fresagem [RPM]	800 1000	800 1000	800 1000	
<b>Fresas</b>	Lança #2,2 <b>5088</b>	Helicoidal #2,2 <b>5073</b>	Cônica #2,8 <b>5084</b>	Countersink SP/RP <b>5047</b>

Osso tipo I, II, III y IV

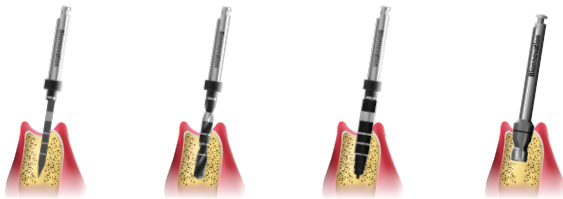

<b>Implantes</b> Diâmetro	3,3	3,3	3,3	Opcional
------------------------------	-----	-----	-----	----------

